


Borderless Lab 365 – 預約系統


1. 前往 <https://stem-ap.polyu.edu.hk/remotelab/>
2. 點擊“立即預約”



3. 選取你的實驗及點擊下方“立即預約” 按鈕

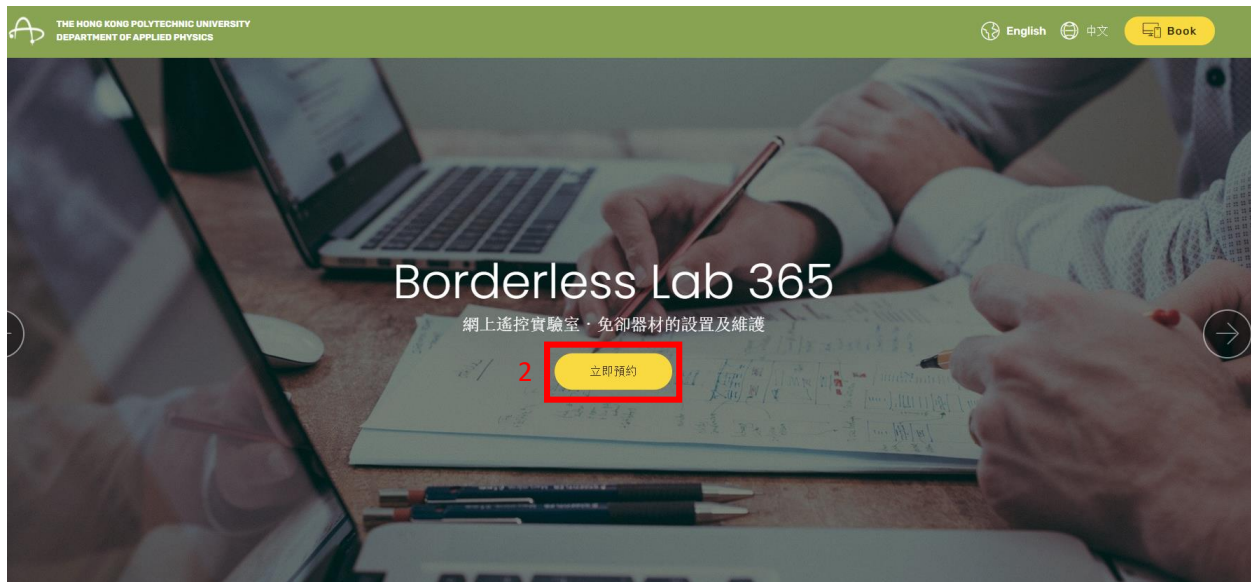

4. 輸入你的登入名稱及密碼

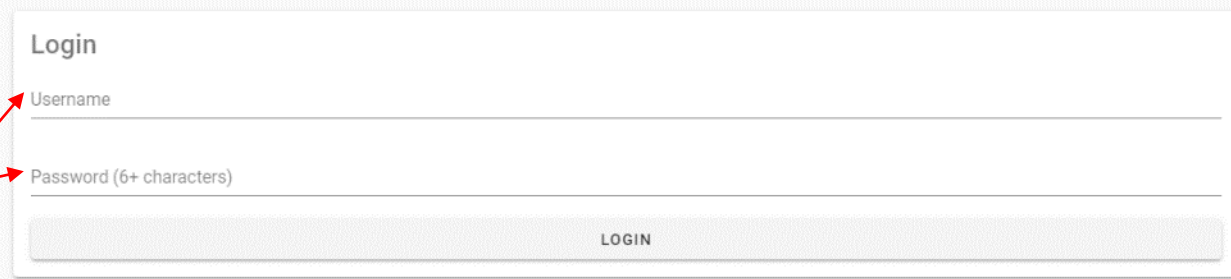

Login

Username

Password (6+ characters)


LOGIN

5. 選取需要預約的時間然後點擊“BOOK”




From	To	Room	Action
15:00	15:55	Room 1	5 BOOK
16:00	16:55	Room 1	BOOK
17:00	17:55	Room 1	BOOK
18:00	18:55	Room 1	BOOK
19:00	19:55	Room 1	BOOK
20:00	20:55	Room 1	BOOK
21:00	21:55	Room 1	BOOK
22:00	22:55	Room 1	BOOK
23:00	23:55	Room 1	BOOK

6. 在預約時間登入系統後，點擊“ENTER LAB”開始實驗

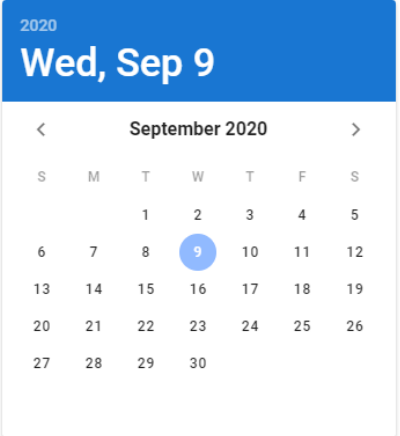

From	To	Room	Action
15:00	15:55	Room 1	6 ENTER LAB
16:00	16:55	Room 1	BOOK
17:00	17:55	Room 1	BOOK
18:00	18:55	Room 1	BOOK
19:00	19:55	Room 1	BOOK
20:00	20:55	Room 1	BOOK
21:00	21:55	Room 1	BOOK
22:00	22:55	Room 1	BOOK
23:00	23:55	Room 1	BOOK

7. 根據實驗手冊指示進行實驗，完成後點擊“LOGOUT”離開

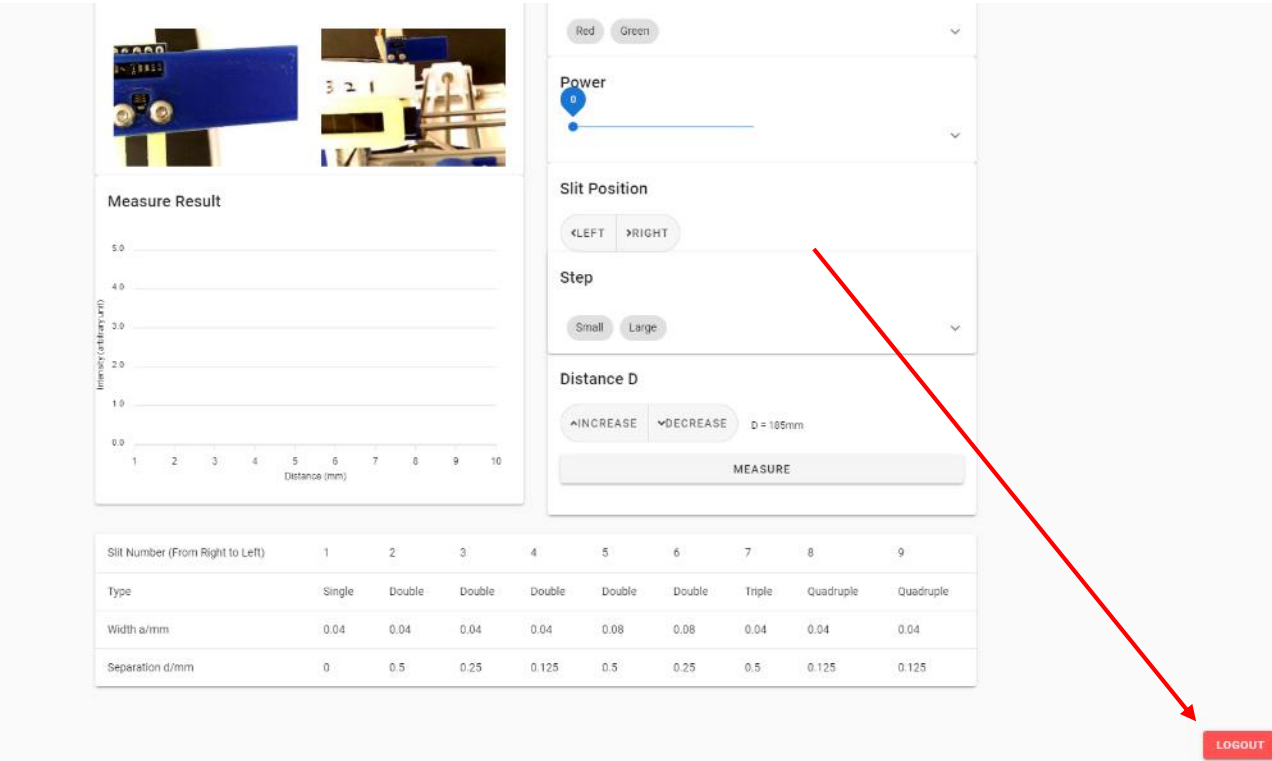

The screenshot displays the experimental control interface for Borderless Lab 365. It includes a 'Measure Result' graph, a table of slit parameters, and various control buttons. A red arrow indicates the path from the 'MEASURE' button to the 'LOGOUT' button.

Measure Result

Intensity (arbitrary unit)

Distance (mm)

Slit Number (From Right to Left)	1	2	3	4	5	6	7	8	9
Type	Single	Double	Double	Double	Double	Double	Triple	Quadruple	Quadruple
Width a/mm	0.04	0.04	0.04	0.04	0.08	0.08	0.04	0.04	0.04
Separation d/mm	0	0.5	0.25	0.125	0.5	0.25	0.5	0.125	0.125

Control Panel:

- Red Green
- Power
- Slit Position: <LEFT >RIGHT
- Step: Small Large
- Distance D: INCREASE DECREASE D = 105mm
- MEASURE
- LOGOUT

8. 如有任何疑問，歡迎電郵至 klaiwong@polyu.edu.hk 與黃小姐聯絡